



Variable-Gain Photoreceiver - Fast Optical Power Meter

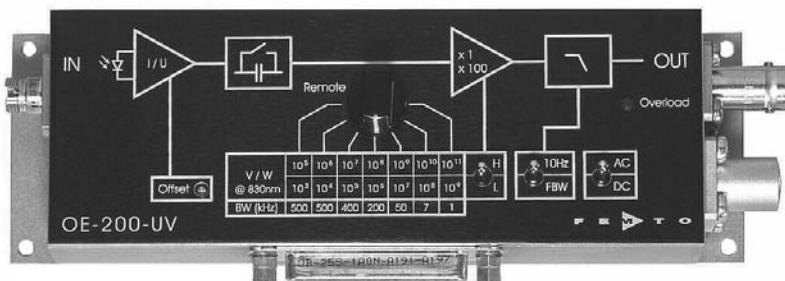

Features	<ul style="list-style-type: none"> Si Detector, 1.1 x 1.1 mm² Active Area, FC, ST or SMA Fiber Receptacle Spectral Range 190 - 1000 nm, UV-Enhanced Conversion Gain Switchable from 1 x 10³ to 1 x 10¹¹ V/W Calibrated at 830 nm, Traceable to NIST Standards Bandwidth up to 500 kHz Local and Remote Control
Applications	<ul style="list-style-type: none"> Fast Optical Power Meter Spectroscopy General-Purpose Opto-Electrical Measurements Optical Receiver for Use with Lock-In Amplifiers

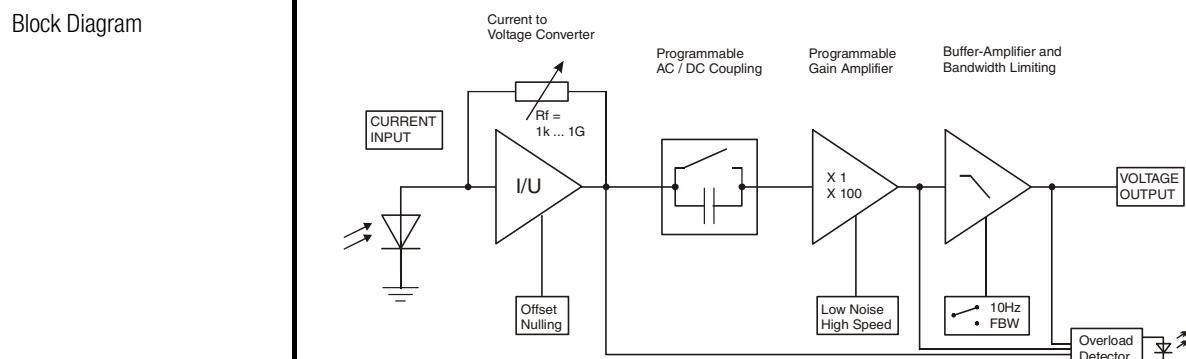


BS01-0149059-

Variable-Gain Photoreceiver - Fast Optical Power Meter

Specifications	<i>Test Conditions</i>	$V_s = \pm 15 V, T_a = 25^\circ C$								
Gain	Conversion Gain Gain Accuracy Gain Drift				$1 \times 10^3 \dots 1 \times 10^{11} V/W$ (@ 830 nm) $\pm 5\%$ electro-optical ($P_{opt} \leq 1 \text{ mW}$), traceable to NIST see table below					
Frequency Response	Lower Cut-Off Frequency Upper Cut-Off Frequency Gain Flatness Gain Accuracy				DC / 1 Hz, switchable up to 500 kHz (See Table), switchable to 10 Hz $\pm 0.1 \text{ dB}$ $\pm 1\%$ electrical, between settings					
Input	NEP Max. cw-Saturation Power Dark Current Compensation				see table see table $\pm 2 \text{ nW}$, adjustable by offset trimpot and external control voltage					
Detector	Detector Active Area Spectral Response Sensitivity Dark Current				Si photodiode in FC, ST or SMA fiber receptacle $1.1 \times 1.1 \text{ mm}^2$ 190 – 1000 nm, UV-enhanced 0.3 A/W (@ 830 nm) 0.1 A/W (@ 200 nm) 2 pA typ.					
Performance Depending on Gain Setting	Gain Setting (Low Noise) (V/W)	10^3	10^4	10^5	10^6	10^7	10^8	10^9		
	Upper Cut-Off Frequency (- 3 dB) Rise / Fall Time (10% - 90%) NEP ($/\sqrt{\text{Hz}}$, @100 Hz) Offset Current Drift ($^\circ\text{C}$) Gain Drift ($^\circ\text{C}$) cw-Saturation Power	500 kHz 700 ns	500 kHz 700 ns	400 kHz 900 ns	200 kHz 1.8 μs	45 kHz 8 μs	7 kHz 50 μs	1.2 kHz 300 μs		
	Gain Setting (High Speed) (V/W)	10^5	10^6	10^7	10^8	10^9	10^{10}	10^{11}		
	Upper Cut-Off Frequency (- 3 dB) Rise / Fall Time (10% - 90%) Min. NEP ($/\sqrt{\text{Hz}}$, @100 Hz) Offset Current Drift ($^\circ\text{C}$) Gain Drift ($^\circ\text{C}$) cw-Saturation Power	500 kHz 700 ns	500 kHz 700 ns	400 kHz 900 ns	200 kHz 1.8 μs	45 kHz 8 μs	7 kHz 50 μs	1.2 kHz 300 μs		
Output	Output Voltage Output Impedance Max. Output Current				$\pm 10 \text{ V}$ (@ $> 10 \text{ k}\Omega$ load) $50 \text{ }\Omega$ (terminate with $> 10 \text{ k}\Omega$ load for best performance) $\pm 30 \text{ mA}$					
Indicator LED	Function	Overload								
Digital Control	Control Input Voltage Range Control Input Current Overload Output	Low: - 0.8 ... + 1.2 V, High: 2.3 ... + 12 V 0 mA @ 0 V, 1.5 mA @ + 5 V, 4.5 mA @ + 12 V non active: 0 V, max. -1 mA, active: 5.1 V, max. 7 mA								
Ext. Offset Control	Control Voltage Range Offset Control Input Impedance	$\pm 10 \text{ V}$ $20 \text{ k}\Omega$								

Variable-Gain Photoreceiver - Fast Optical Power Meter

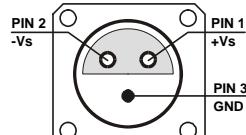
Specifications (continued)		
Power Supply	Supply Voltage Supply Current Stabilized Power Supply Output	± 15 V + 150 / -100 mA (depends on operating conditions, recommended power supply capability minimum 250 mA) ± 12 V, max. 150 mA, + 5V, max. 50 mA
Case	Weight Material	320 g (0.74 lbs) AlMg4.5Mn, nickel-plated
Temperature Range	Storage Temperature Operating Temperature	-40 ... +80 °C 0 ... +60 °C
Absolute Maximum Ratings	Max. cw-Power (averaged) Control Input Voltage Power Supply Voltage	20 mW - 5 V / + 16 V ± 22 V
Connectors	Input Output Power Supply Control Port	optical, FC, ST or SMA fiber receptacle BNC LEMO series 1S, 3-pin fixed socket Pin 1: + 15V Pin 2: - 15V Pin 3: GND  Sub-D 25-pin, female, qual. class 2 Pin 1: +12V (stabilized power supply output) Pin 2: -12V (stabilized power supply output) Pin 3: AGND (analog ground) Pin 4: +5V (stabilized power supply output) Pin 5: digital output: High = overload Pin 6: signal output (connected to BNC) Pin 7: NC Pin 8: input offset control voltage Pin 9: DGND (ground for digital control pins 10 - 14) Pin 10: digital control input: gain, LSB Pin 11: digital control input: gain Pin 12: digital control input: gain, MSB Pin 13: digital control input: AC/DC Pin 14: digital control input: high speed / low noise Pin 15 - 25: NC
Available Models	OE-200-UV-FC OE-200-UV-ST OE-200-UV-SMA	FC receptacle, calibrated at 830 nm ST receptacle, calibrated at 830 nm SMA receptacle, calibrated at 830 nm

Variable-Gain Photoreceiver - Fast Optical Power Meter

Remote Control Operation	General	<p>Remote control input bits are opto-isolated and connected by a logical OR function to the local switch settings. For remote control set the corresponding local switches to "Remote", "AC" and "H" and select the desired setting via a bit-code at the corresponding digital inputs.</p> <p>Mixed operation, e.g. local AC/DC setting and remote controlled gain setting, is also possible.</p> <p>The switch setting "FBW / 10 Hz" of the lowpass signal filter is not remote controllable.</p>				
	Gain Setting	Low Noise Gain (V/W) Pin 14=High	High Speed Gain (V/W) Pin 14=Low	Pin 10 LSB	Pin 11	Pin 12 MSB
		10^3	10^5	Low	Low	Low
		10^4	10^6	High	Low	Low
Spectral Response	AC/DC Setting	Coupling	Pin 13			
		AC	Low			
		DC	High			

Variable-Gain Photoreceiver - Fast Optical Power Meter

Dimensions

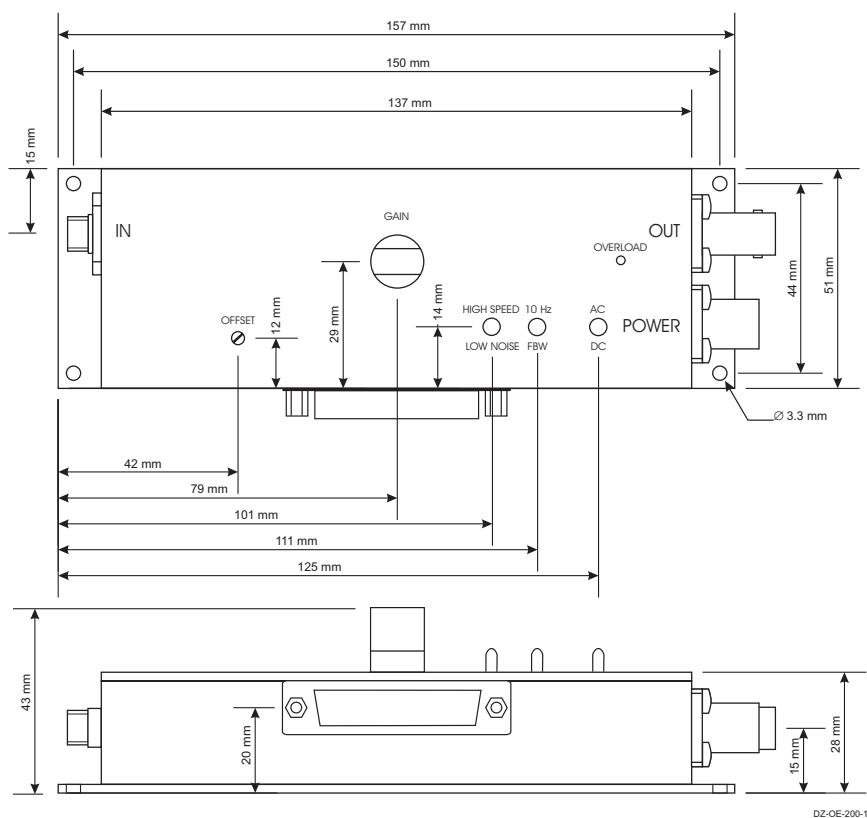

FEMTO Messtechnik GmbH
Paul-Lincke-Ufer 34
D-10999 Berlin • Germany
Tel.: +49 (0)30 – 4 46 93 86
Fax: +49 (0)30 – 4 46 93 88
e-mail: info@femto.de
<http://www.femto.de>

Specifications are subject to change without notice. Information furnished herin is believed to be accurate and reliable. However, no responsibility is assumed by FEMTO Messtechnik GmbH for its use, nor for any infringement of patents or other rights granted by implication or otherwise under any patent rights of FEMTO Messtechnik GmbH. Product names mentioned may also be trademarks used here for identification purposes only.

© by FEMTO Messtechnik GmbH

Printed in Germany